

# BARI10A


## BARIATRIC BED



### Staff Control Operation:

- 1 Push button to raise back
- 2 Push button to lower back
- 3 Push button to raise knee
- 4 Push button to lower knee
- 5 Push button to raise bed
- 6 Push button to lower bed
- 7 Push button for comfort chair position
- 8 To disengage chair mode, press and hold button #8 until bed is completely flat
- 9 Push button for Trendelenburg
- 10 Push button for Reverse Trendelenburg
- 11 Push and hold to change lockout status in combination with another function button
- 12 Push button for staff control backlighting

**Note:** Use hand pendent for all controls when footboard controls are not present.

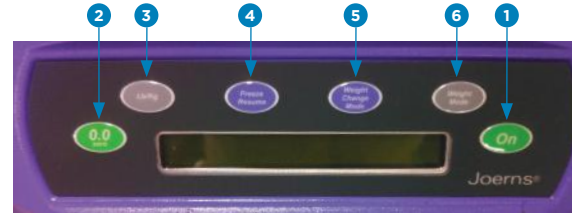

### Scale System Function:

- 1 **Power:** Press ON. Patient's weight will be shown on screen.
- 2 **Zero the Scale:** With only linens on bed, press 0.0 ZERO button. "PLEASE WAIT-HANDS OFF" will display. (Note: Do not touch the weight-sensitive portion of integrated scale system.) In a few seconds, display will read "WEIGHT = 0.0 LB".
- 3 **Changing between lbs/kgs:** Press LB/KG to easily switch between the two measurements.
- 4 **Freeze/Resume:** To add/remove equipment without affecting the patient's weight reading, press the Freeze/Resume mode button and display will read "FREEZE = (patient's weight)." Make all necessary adjustments to bed or patient (add/remove any pillows, blankets, equipment, etc). When finished, press Freeze/Resume button again to resume normal patient weighing.
- 5 **Weight Change Mode:** To monitor patient's weight relative to zero, press the Weight Change Mode button, which records weight variations.
- 6 **Weight Mode:** Press the Weight Mode button to view patient's weight in either lbs or kgs.

# BARI10A

## BARIATRIC BED



**Length Extension:** Slide the frame extender out by pulling the latch knob on the side of the bed frame. Next, slide the deck extender by pulling the latch knob underneath the deck and sliding the extender out until the latch pin clicks into its detent.



**Width Extension:** Pull purple latch handle along the bottom edge of each section and pull.

### Emergency Release:

To lower bed to a flat position, use pendant and press down arrow for each option.

In case of power or battery backup failure, turn red emergency release lever up and hold in place until platform lowers to desired position; release handle to stop descent.



### Bed Reset (bed alarms or freezes):

1. Buttons may turn an orange hue.
2. Press and hold buttons 5 and 6 at the exact same time, the bed should start dinging. If it does not, you must retry this step.
3. Listen for 8-9 dings.
4. Hold buttons until dinging stops.
5. Release buttons.
6. Press button 6 to lower the bed to the lowest position.
7. Once in lowest position the bed will be completely reset.



To learn more about Agiliti solutions go to [www.agilitihealth.com](http://www.agilitihealth.com). For further questions, refer to Bari10A User Manual or contact Agiliti at **1.800.847.7368**.