


NP Adjust™

Self-Adjusting Support Surface
with CuPro™ foam



CuPro™ foam

NP Adjust features self-adjusting foam-filled air cells which respond to patient weight and movement—promoting immersion, pressure redistribution and comfort. Four-way stretch top cover fabric and a layer of CuPro copper- and gel-infused foam help enhance microclimate by dissipating heat and moisture. Highly chemically resistant performance fabrics, RF-welded seams and exclusive Core Shield™ layer protect the internal components from contamination.

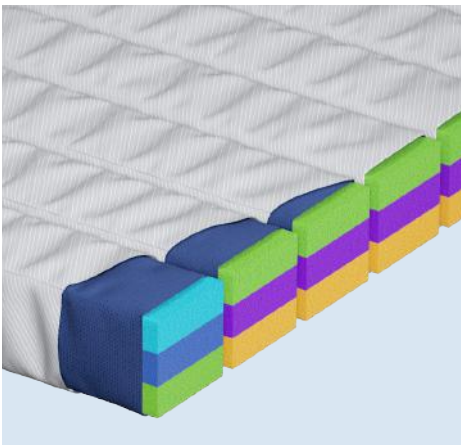

Self-adjusting cells contain layered and zoned foam for pressure redistribution



4-way stretch, high-performance breathable top cover promotes immersion



Core Shield™ minimizes risk of fluid ingress, extending surface life

KEY SPECS

Max Patient Weight 500 lb.
Width 35"
Length 82"
Height 7"

CONSTRUCTION DETAILS

Foam-Filled Air Cells Standard
Edge Support Standard
RF-Welded Seams Standard

OPTIONS

Behavioral Health NP Adjust BH™
Seamless Expandable NP Adjust EX™
Custom Sizes By Request

LEARN MORE

Ask your Agiliti representative or call 800-814-9389

