


ELITE ACUTE CARE BED WITH POWER ASSIST DRIVE

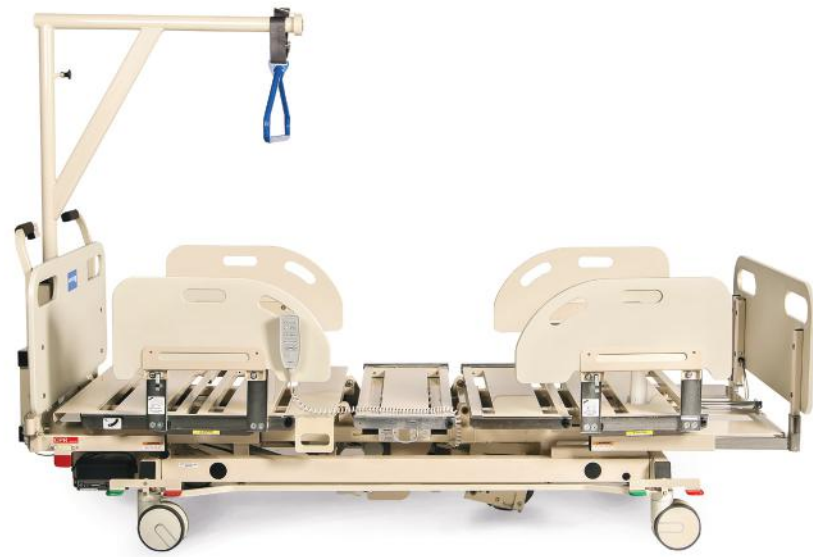

PATIENT WEIGHT CAPACITY: 1,000 pounds

FOWLER BACK REST ANGLE: 0-65 degrees

FOLEY BAG HOLDER: One attachment each side of bed

RESTRAINT ATTACHMENT: One attaching point, each side of bed

POWER ASSIST DRIVE: Engage P.A.D. by lowering bed deck to lowest position. Disengage by raising bed deck

VARIABLE SPEED CONTROL: Switch embedded in right side head end. Push handle for forward/reverse Power Assist Drive functions

TRENDELEBURG/REVERSE TRENDELEBURG: 9 degrees

KNEE FLEX/GATCH: 0-20 degrees

DECK LOW POSITION: 17"

DECK HIGH POSITION: 26"

ADJUSTABLE DECK LENGTH: 84"-92"

ADJUSTABLE DECK WIDTH: 39"-48"-54"

SIDE RAILS (4): 17" high above deck

CPR FOWLER BACKREST RELEASE: Head end, left and right

FOUR WHEEL LOCK AND STEER: Engage wheel locks or steer at any caster wheel

INTEGRATED PUSH HANDLES: Two, at head end

PATIENT HELPER/TRAPEZE: Swings 90 degrees from center

CONTROL PENDANT: Illuminated, five-function, self-storing

CAREGIVER CONTROL: Controls all functions and locks out pendant functions

BED POSITIONING/ARTICULATION: High/low, fowler backrest, Trendelenburg/Reverse Trendelenburg, knee flex/leg lift, cardio chair

ANGLE INDICATOR: Two at head end, one left, one right

UNDER BED ILLUMINATION: LED illumination

OUT OF BED EXIT ALARM: Sounds at bed with option to sound at nurse call station using DB 37 cable or 1/4" phone connector

INTEGRAL SCALE WEIGHING SYSTEM:

- Displays accurate weight in any articulated position including Trendelenburg and Reverse Trendelenburg and in any caster orientation
- One convenient control box contains both the scale controls and the caregiver key pad
- Scale display rotates - HIPAA compliant
- Powered by bed electrical system and battery backup. Does not require batteries and will function when disconnected from A/C power
- Weight change mode
- Dual weight display
- Freeze mode
- Backlit display
- Accurate to within 0.2% or 2 pounds, whichever is greater

PATENT PENDING
DESIGNED and MANUFACTURED IN USA

WHY BUY GENDRON:

- 25 years experience serving the bariatric patient
- Designed and manufactured in the USA
- Complete bariatric product portfolio



800.537.2521 • sales@gendroninc.com

www.gendroninc.com

GENDRON

TOTAL
SOLUTIONS
FOR
BARIATRIC
PATIENT
CARE

MOBILITY
SAFETY
COMFORT
DIGNITY

The Elite Acute Care Bed with Power Assist provides your bariatric patients with all of the above in an easy-to-use design, engineered to allow you to give your patients the best possible care. Now, with the addition of Power Assist Drive, the Elite bed provides maximum mobility and ease of use for staff and transport team. Easily navigate and maneuver with the Power Assist Drive.



ELITE BARIATRIC ACUTE CARE BED
FEATURING POWER ASSIST DRIVE

Elite Acute Care Bed with Power Assist Drive Features

VARIABLE SPEED SWITCH, FORWARD AND REVERSE

ON/OFF POWER SWITCH FOR POWER ASSIST DRIVE

OPTIONAL LOW PROFILE PATIENT HELPER
Transfer, lift and position patients with ease, swings 90 degrees from center for accessibility to the patient and eliminates the need to remove from the bed when transferring. The low profile provides unobstructed patient access when the headboard is removed

INTEGRATED BED SCALE SYSTEM
Ensure patient privacy - scale display window rotates to block the view of the information shown on the scale display

CAREGIVER CONTROL
Located at foot end of bed. Any or all hand pendant functions can be locked-out at the caregiver control to restrict patient access

FEATURES INCLUDE:

BED POSITIONING/ ARTICULATION
Provide proper clinical positioning from Trendelenburg and Reverse Trendelenburg to knee flex and cardio chair. Easy-to-use button controls; one button cardio chair positioning and one button return to level positioning

OUT OF BED EXIT ALARM
Alerts you when patient attempts to exit bed. Programmable alarm sounds at the bed, with an option to connect to the nurse station

ADJUSTABLE DECK LENGTH
Ensure the comfort of your tall patients - one touch extension from 84" to 92" long

QUICK ELECTRICAL DISCONNECT FOR HEADBOARD REMOVAL

NURSE CALL CONNECTIONS

CPR FOWLER BACKREST RELEASE
Respond faster to patient needs: lower backrest to flat within six seconds for CPR

FOWLER BACKREST
Minimum 65 degree angle in raised position

FOUR WHEEL LOCK AND STEER
Steer safely and precisely - engage wheel lock and steer at any caster wheel

CONTROL PENDANT
Adjust settings confidently even in dark rooms. Illuminated controls; five-function, self-storing at head end, left and right side

WRIST RESTRAINT ATTACHMENT POINT
Located on both sides of the bed

FOLEY BAG HOLDER
Located on both sides of the bed

ADJUSTABLE DECK WIDTH
Transport safely and provide comfort - one touch extension from 39" to 48" to 54" wide

POWER ASSIST DRIVE MECHANISM
Incorporates mechanical electric brake manual release

DUAL HAND HOLDS
On all four side rails for patient safety