

NPT3™

Gel-Infused Foam Support Surface




Technology x3: top cover, MAP Liner™, and CuPro™ copper-infused foam. Shown with optional alternating pump.

NPT3™ is constructed of 12 air chambers that provide three separate zones of support, with head-to-heel coverage for patient comfort. Vented technology creates expedited pressure equalization to maximize patient immersion and provide advanced pressure redistribution for all body types. NPT3 is an Envy® product, which reduces chemicals of concern in the fabric, foam, and fire barrier.

- RF-welded First Point of Contact™ top cover is free of harmful chemicals
- Mattress Asset Protection (MAP) Liner™ safeguards internal components
- CuPro™ copper-infused, thermoresponsive foam keeps the patient cool while reacting to temperature and weight
- Pump available for added 1-in-2 alternation therapy
- Available for most hospital bed frames

sizewise.com 800-814-9389


Sizewise®
Designed to heal

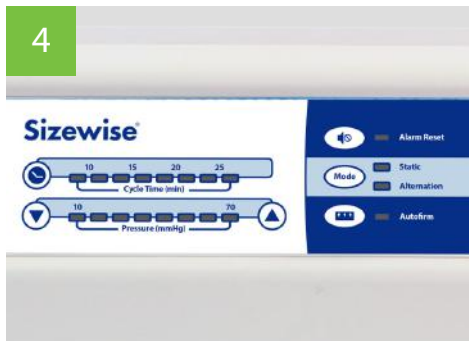

Key Specs

MATTRESS

Capacity	500 lb./1,000 lb.
Mattress Weight	35 lb.
Width	35", 39", or 48"
Length	80-86"
Custom Sizes	Available
Reinforced Evacuation Handles	Standard
Adapt™ Heel Zone	Optional
Alternating Air Pump	Optional

FABRIC OPTIONS

Extreme Polycarb 4-Way Stretch Top	Standard
Non-Skid Polyurethane Bottom Cover	Standard
Halogen-Free Fire Barrier	Standard
MAP Liner™	Standard

FABRIC CHARACTERISTICS



4-Way
Stretch



Chemical
Resistant



Green

- 1 — CuPro™ copper-infused foam keeps the body cool while reacting to temperature and weight
- 2 — Mattress Asset Protection (MAP) Liner™ safeguards internal components and reduces shear
- 3 — Extreme Polycarbonate 4-Way Stretch Top Cover is highly durable against cleaning protocols
- 4 — Optional pump adds 1-in-2 alternation therapy
- 5 — Available to fit most facility-owned bed frames



Learn More:

Visit sizewise.com/NPT3

sizewise.com 800-814-9389


SizeWise®
Designed to heal