

Pediatric Pulse™

Low Air Loss Pulsation Therapy




Available in multiple sizes to
fit cribs and Med-Surg beds

Pediatric Pulse™ is a latex-free mattress replacement system designed for pressure injury treatment and prevention in infants, children, and adolescents. The top cover is vapor-permeable, low-friction, and anti-shear. The three-zoned surface features vented air cells to provide pulsation or static therapy while managing the patient's microclimate. Mattress comes in multiple standard and custom sizes to fit everything from cribs to Med-Surg beds (details on reverse).

- Adjustable pulsation times from three to 20 minutes
- Vapor-permeable, low-friction, anti-shear top cover
- Child-friendly control unit keypad design
- Autofirm feature for safe and comfortable transport
- Available in standard and custom sizes

sizewise.com 800-814-9389


Sizewise®
Designed to heal

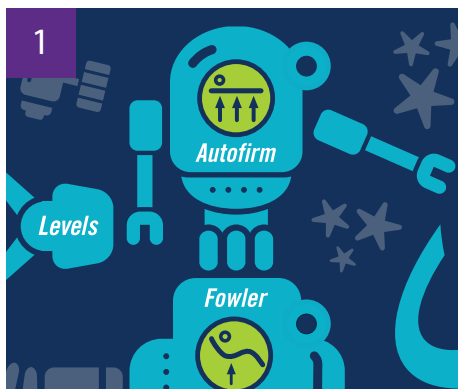

Key Specs

MATTRESS (CRIB)

Min. Patient Weight 1.6 lbs.
 Max. Patient Weight 100 lbs.
 Sizes Standard and Custom

MATTRESS (MED-SURG)

Min. Patient Weight 35 lbs.
 Max. Patient Weight 200 lbs.
 Width 35"
 Length 80"
 Height 8"

CONTROL UNIT

Zoned Low Air Loss Yes
 Pulsating Low Air Loss 3-20 Min.
 Autofirm Yes
 Fowler Yes
 Alarms Power Failure and Low Pressure
 Pressure Sensor Technology Yes
 CPR Release Yes
 Unit Weight 11 lbs.

FABRIC OPTIONS (CRIB)

Nylon 70 Denier Polyurethane-Coated ... Standard
 Extreme Polycarb 4-Way Stretch Optional
 Non-Skid Polyurethane Bottom Cover ... Standard

FABRIC OPTIONS (MED-SURG)

Nylon 70 Denier Polyurethane-Coated ... Standard
 4-Way Stretch Optional
 PVC Laminate Bottom Cover Standard

FABRIC CHARACTERISTICS




Powered
Surface



Cool
Comfort



Low
Friction

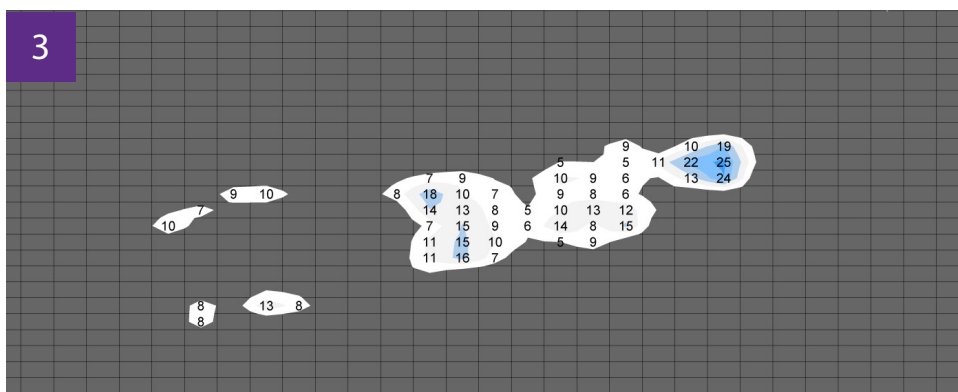

1 — Autofirm feature for safe and comfortable transport

2 — Easy-to-use controls in a child-friendly motif

3 — Effective pressure injury prevention and treatment for patients as small as 1.6 lbs. (male patient, 3 yo, 34.8 lbs. shown)

4 — Available in crib and Med-Surg sizes to fit a variety of patient needs

sizewise.com 800-814-9389


Sizewise[®]
 Designed to heal