


## Quartet™ Mattress



- ▶ Three mattress models available with the following therapies:
  - Auto immersion
  - Low air-loss
  - Pulsation
  - Alternating pressure, turn assist or continuous lateral rotation
- ▶ Active sensors to monitor and automatically adjust immersion settings to match patient weight and position
- ▶ Internal dual battery for power outages and patient transfers
- ▶ Standard and bariatric mattress options

## The Quartet™ Mattress Series: Customizing the care your patients need

The OnCare™ Quartet™ solves the challenge of managing high-acuity, complex wounds in one easy-to-use, comprehensive package. Engineered with intelligent auto immersion software that continuously monitors and adjusts to a patient's changing positions, the Quartet also offers low air-loss, pulsation and your choice of alternating pressure, turn assist or continuous lateral rotation. These features not only accelerate the recovery process, but reduce the burden on clinicians to continuously monitor patient movements and mattress settings, allowing for more time to spend on bedside care.

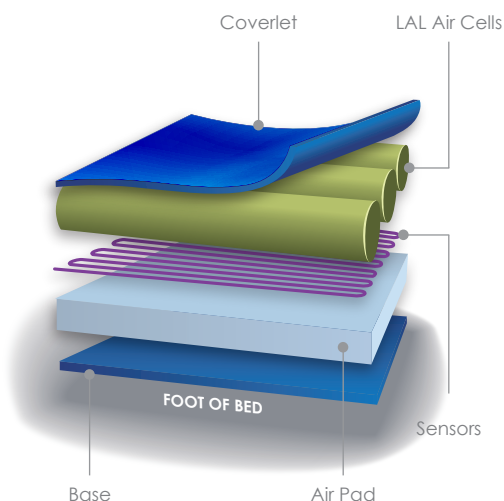



# Quartet™ Mattress

## Features and Specifications

### Clinical Applications:

- Stage I through IV pressure ulcers
- High-risk patients requiring aggressive wound prevention
- Excessive maceration
- Progressive motor decline
- Acute surgical wounds

Therapy Features	Quartet Mattress Models		
	Auto Immersion	Turn Assist	Clear Turn
Auto Immersion	✓	✓	✓
Low Air-Loss	✓	✓	✓
Pulsation	✓	✓	✓
Alternating Pressure	✓		
Heel Zone Regulator	✓		
Turn Assist		✓	
Continuous Lateral Rotation			✓

### Standard Features for All Models

- Max Inflate
- Internal Dual Batteries
- Lock Out
- Alarms
- CPR Deflate
- Five Zone Mattress
- Removable Air Cells Allow for Configurable Surface
- Fowler Boost
- Custom Comfort Settings
- Bottom Out Protection
- Breathable Non-Shear Nylon Top Sheet
- Cable Management System

### Specifications:

	Size	Type	Weight
Quartet™ Standard Replacement Mattress	80"L x 36"W x 10"H	Flat Deck and Recessed Deck	Capacity: 500 lbs.
	84"L x 36"W x 10"H		
Quartet™ Bariatric Replacement Mattress	80"L x 48"W x 10"H	Flat Deck	Capacity: 1000 lbs.
	84"L x 48"W x 10"H		

*Note: Additional size options are available for custom order*

CORPORATE OFFICE 6625 West 78<sup>th</sup> Street, Suite 300 Minneapolis, MN 55439  
Toll Free: 800.847.7368 Phone: 952.893.3200 Fax: 952.893.0704 www.uhs.com



OnCare™ Quartet™ is a proud offering of Universal Hospital Services, Inc.  
Patented Product US #8102270, US #7454809, US #8220090 and UHS #7171711.  
©2016 Universal Hospital Services, Inc. All rights reserved. QUARTET.0316

