


# ROTATION LOW AIR LOSS

to promote  
pulmonary hygiene

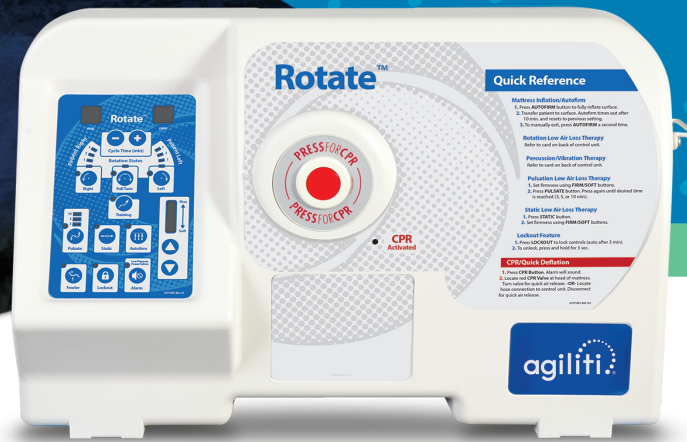

## Rotate™/Big Turn™

Rotate and its bariatric counterpart, Big Turn, provide continuous lateral rotation therapy (CLRT) and immersion for pressure redistribution. The surfaces' low air loss capability and high-tech fabric combine to help control patient microclimate and both surfaces feature a pulsation mode, which stimulates capillary blood flow to the skin. Rotate and Big Turn are sized to fit most Med-Surg and ICU bed frames.



Rotation, immersion,  
and pulsation modes



Optional percussion/  
vibration module



Ideal for immobile  
patients at risk for  
pressure injury



Breathable, low-friction,  
anti-shear top  
cover fabric



**LEARN MORE**


[agilithhealth.com](http://agilithhealth.com)

800-814-9389




# Easy-to-Use Interface

includes instructions for use and QR code for training/support resources

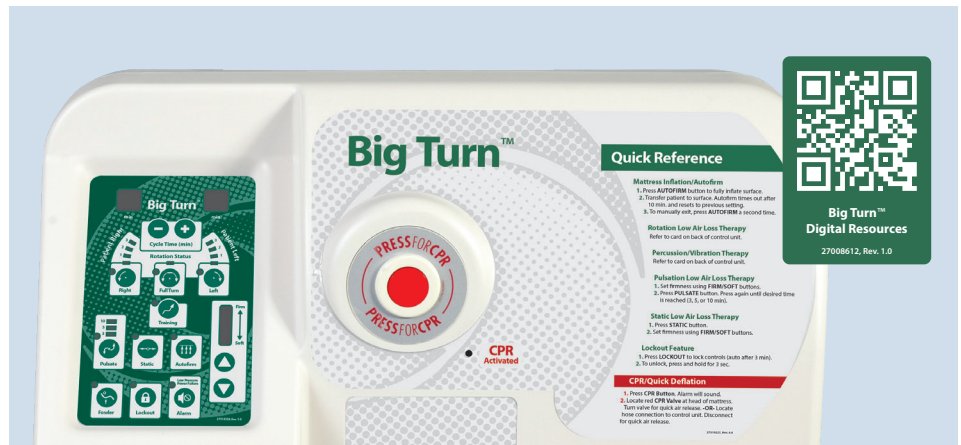

Full-body continuous lateral rotation therapy (CLRT), including left/right/full turn and turn assist

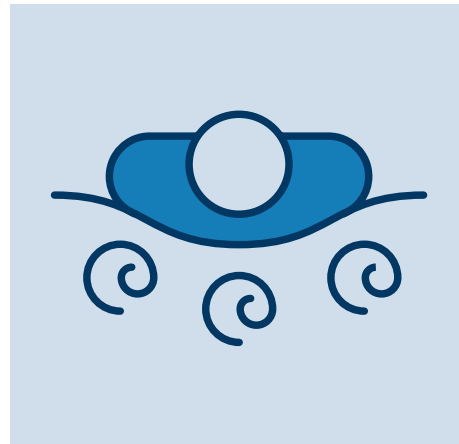

Low air loss technology and breathable top cover help control patient microclimate



Optional percussion/vibration module helps loosen secretions in the lungs

## KEY SPECS

### MATTRESS

Safe Working Load . . . . .600 lb./1,000 lb.  
 Max Patient Weight . . . . . 550 lb./950 lb.  
 Max Turning Weight . . . . . 450 lb./750 lb.  
 Width . . . . . 35", 39", or 48"  
 Length . . . . . 80" or 82"  
 Height . . . . . 10"  
 Weight . . . . . 19-37 lb.  
 Sizes . . . . . Standard and Bariatric

### CONTROL UNIT

Rotation Mode . . . . . 1-30 min.  
 Pulsation Mode . . . . . 3, 5, or 10 min.  
 Immersion Mode (Static) . . . . . Yes  
 Training/Acclimation Mode . . . . . Yes  
 Low Air Loss . . . . . Yes  
 Alarms . . . . . Power Failure/Low Pressure  
 Control Unit Type . . . . . Blower

CPR Release . . . . . Yes  
 HEPA Filter . . . . . Yes  
 Lockout . . . . . Yes  
 Percussion/Vibration . . . . . Optional  
 Pressure Sensor Technology . . . . . Yes  
 Turn Assist . . . . . Yes  
 Weight . . . . . 16 lb.

### PERFORMANCE FABRICS

Polycarbonate 4-Way Stretch . . . . .  
 . . . . . Standard Top Cover  
 Nylon Taffeta . . . . .  
 . . . . . Optional Top Cover  
 Non-Skid Polyurethane . . . . .  
 . . . . . Standard Bottom Cover

### FABRIC CHARACTERISTICS



Breathable



Low Friction



Advanced Chemical Resistance