

Secure II® Med/Surg Bed

Available with Stryker's patented Chaperone® bed exit system and Zone Control® technology

Single-button, in-bed scale provides accurate, repeatable weights with the patient in any position

Open architecture allows the flexibility to choose any sleep surface

Retractable bed frame helps keep patients close to their belongings regardless of bed articulation

Fall prevention educational program and CEU course available

iBed Awareness™ option monitors brake, siderail position, bed exit and low bed height, and provides proactive alert if parameters change



The Standard in Patient Care

Responsive engineering over time has resulted in a med/surg bed that offers the highest level of security, intuitive operation and flexible functionality. Fully-customizable through a host of features designed around patient and caregiver satisfaction, the Secure II represents the benchmark in patient handling.


The combination of Stryker's Chaperone bed exit system and Fall Prevention Program, coupled with the feature of bed retractability, form the most comprehensive fall prevention solution available.



Optional Zoom® drive system minimizes the effort needed for patient transport while encouraging proper ergonomic posture.



Intermediate siderail position and low bed height assist with patient entry and exit.

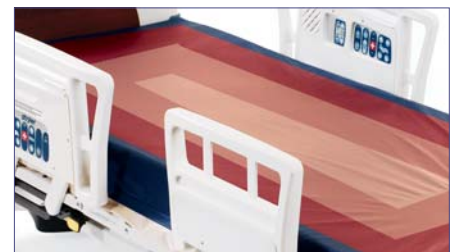

Chaperone accurately senses body positioning by constantly tracking the patient's center of gravity.

Secure II[®]

Med/Surg Bed

stryker[®]



Medical

3800 E. Centre Ave.
Portage, MI 49002
t: 269 329 2100 f: 269 329 2311
toll free: 800 787 9537

www.stryker.com

Standard Features

- 6" (15 cm) casters
- Dual pedestal design
- CPR release
- In-rail fowler/knee gatch control
- Electronic function lock-out control
- Four-wheel, steel-ring brake system with centrally-located activation pedals
- Backrest angle indicator
- Locking-caster steering
- Four traction sockets
- Removable headboard with CPR capability
- Retractable litter top
- Four foley bag hooks
- Four IV receptacles
- Patient restraint locations
- Trend./reverse Trend. controls

EX Model: All Standard Features Plus

- Chaperone[®] center-of-gravity bed exit system
- Scale with isolated foley bag hooks and Trend. display
- Stryker pendant port

ZX Model: EX Features Plus

- Zoom[®] Motorized Drive System

Specifications

• Model Number		
Domestic	3002 (115 VAC)	
International	3221 (230 VAC)	
• Overall Length	93"	(236 cm)
• Overall Width		
Siderails Up	42.5"	(108 cm)
Siderails Down	40"	(102 cm)
• Weight Capacity	500 lbs	(228 kg)
• Height Range (to litter top)		
High	30"	(76 cm)
Low (standard with 6" casters)	16"	(41 cm)
Low (ZX option)	19.75"	(50 cm)
• Litter Positioning		
Backrest	0° - 60°	
Knee Gatch	0° - 40°	
Trend./reverse Trend.	±12° (± 10° with ZX option)	
Retraction	14"	(36 cm)
• Patient Surface	35" x 84"	(89 x 213 cm)
• Caster Diameter		
Standard	6"	(15 cm)
Optional	8"	(20 cm)

Electronics

• Agency Approvals	C-US UL (3002), CE (3221)
• Volts	115 VAC (Domestic) / 230 VAC (International)
• Ampere Rating	7A (Domestic) / 4A (International) 10A (rating for 115 volt outlet option)
• Current Leakage	<100 microamperes
• Frequency	60 HZ (Domestic) / 50-60 HZ (International)
• Hospital Grade Plug	Yes

Stryker reserves the right to change specifications without notice.

Optional Features

- In-bed scale system
- Chaperone center-of-gravity bed exit system
- Chaperone with Zone Control[®]
- Siderail communications—includes nurse call with speaker, TV, radio, volume, room/read lights
- Smart TV—includes closed caption, channel up/down and mute
- Mattresses: management and prevention
- Adapter bracket for patient helper (not available with Zoom)
- Bed extender
- Removable CPR board
- Defibrillator tray
- Pendant control
- Pillow speaker interface
- Pump holder mount on footboard
- Siderail pads
- Transducer mount—siderail or IV pole
- Traction sleeve adapter kits
- Upright oxygen bottle holder
- Wallsaver[™] quick-release connecting kit
- 8" (20 cm) Omni Surface[™] casters with wheel covers
- IV poles
 - Permanent folding two-stage
 - Removable
- 110 volt power outlet (not available with Zoom)
- Prop rod
- iBed Awareness[™] Option

Copyright © 2005 Stryker
LSS-453-06E
PC-SCG-7.5M-1206
Mkt Lit-78-02022005 Rev E