

Adapt Air™


Alternating Pressure/Immersion
Air Support Surface



Adapt Air™ offers pressure injury prevention and treatment through two clinically effective therapy modes: 1-in-2 alternating pressure or powered immersion. In either mode, targeted airflow microclimate management provides a constant flow of air to the support surface's top cover and proprietary AirSpace™ mesh — helping keep patients cool and dry. The quiet, easy-to-use pump provides eight levels of customizable support and can be used on all Adapt line support surfaces.



Performance fabrics and RF-welded seams boost durability and prevent fluid ingress

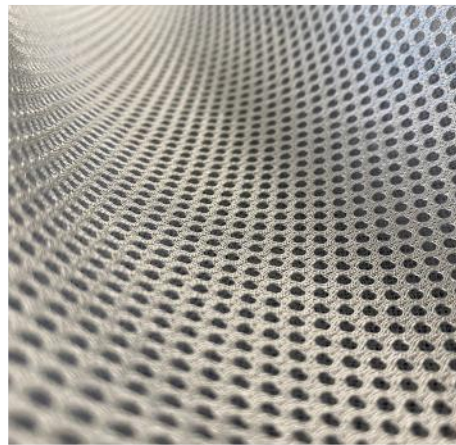

AirSpace™ mesh maintains airflow while reducing friction and shear



Easy-to-use pump allows for therapy selection and comfort adjustment

KEY SPECS (Standard/EX)

Max Patient Weight 550 lb./1,000 lb.
Width 35"/35"-48"
Length 82"/82"-86"
Height: Air Cells/Total 6"/8"

ADAPT PUMP™

Weight 11 lb.
Alternating Cycle Times 5/10/15 min.
Immersion Mode Standard
Microclimate Management Standard


OTHER DETAILS

Alarms Power Failure/Low Pressure
Comfort Settings 8
Seamless Expandable Adapt Air EX™
Transport Cap. Standard

LEARN MORE

Ask your Agiliti representative or call 800-814-9389



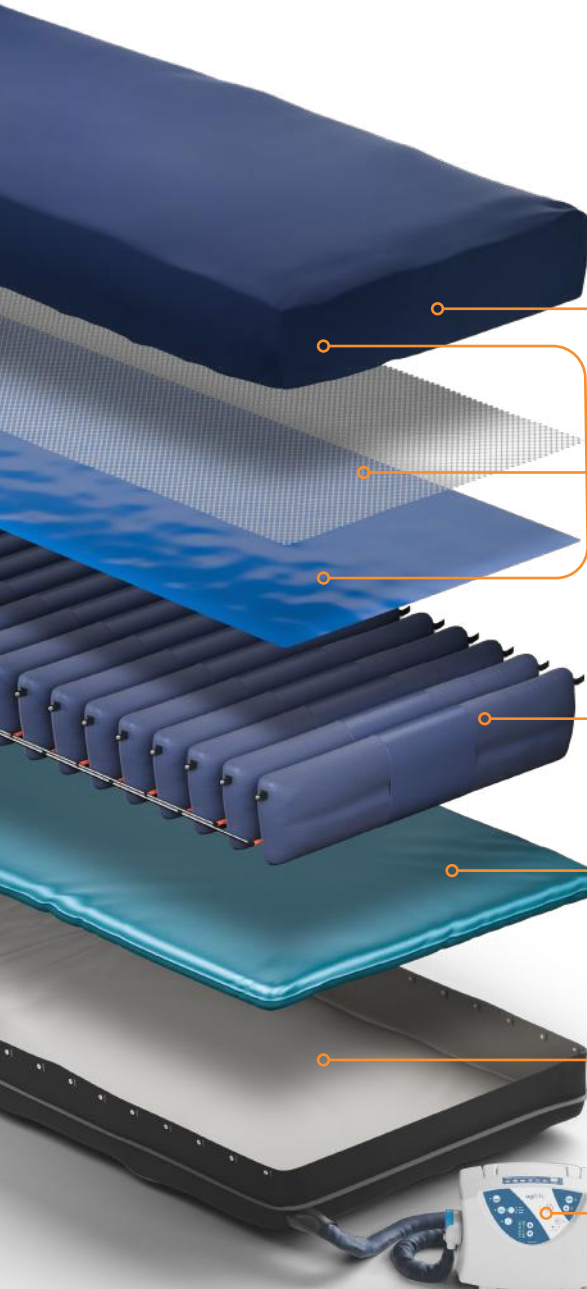

LAYER BY LAYER: ADAPT AIR™

The **4-way stretch top cover** promotes immersion and envelopment; is breathable and chemically resistant, standing up to the toughest cleaning protocols.

Targeted airflow microclimate management combines the top cover fabric, Adapt Pump™ and **AirSpace™ mesh** to ensure non-stop airflow beneath the patient/surface interface — helping keep patients cool and dry.

Twenty 6-inch air cells provide alternating pressure or powered immersion. Unique integrated I-beams maintain vertical structure when weight is applied, leading to deeper immersion potential and lower shear forces.

Though the surface can maintain inflation without power, a **layer of foam** provides additional protection to help prevent bottoming out.

The **bottom cover** is RF-welded from a single panel of non-skid material to minimize the risk of fluid ingress and improve its durability.

The quiet, easy-to-use **Adapt Pump** is compatible across all Adapt line support surfaces.

End-to-End Bed Services


Agility is the only manufacturer of therapeutic support surfaces with end-to-end services that support clinicians and patients, drive measurable outcomes and reduce cost.



Quick access to therapy options, both rental and purchase



Expedite time-to-therapy with onsite management and logistics

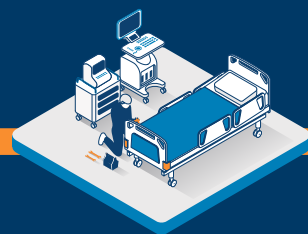

Optimize useful life with comprehensive maintenance and repair

Learn more at agilityhealth.com

PRD-0067C-0524 © 2024 Agility Health, Inc.

