

Adapt Air™


Alternating Pressure/Immersion
Air Support Surface



Adapt Air offers pressure injury prevention and treatment through two clinically effective therapy modes: 1-in-2 alternating pressure or powered immersion. In either mode, Microclimate Assist™ provides a constant flow of air to the support surface's top cover and proprietary AirSpace™ 3D mesh—helping keep patients cool and dry. The quiet, easy-to-use pump provides eight levels of customizable support and can be used on all Adapt Line™ support surfaces.



Performance fabrics and RF-welded seams boost durability and prevent fluid ingress

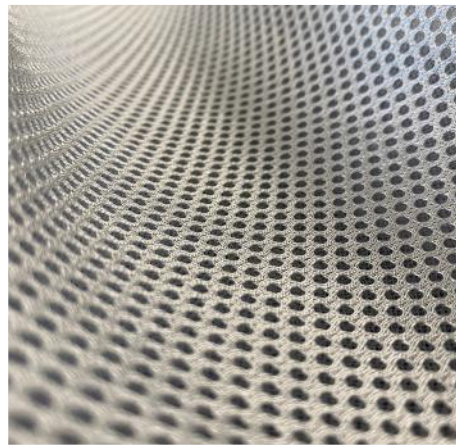

AirSpace 3D mesh maintains airflow while reducing friction and shear



Easy-to-use pump allows for therapy selection and comfort adjustment

KEY SPECS (Standard/EX)

Max Patient Weight 550 lb./1,000 lb.
Width 35"/35"–48"
Length 82"/82"–86"
Height: Air Cells/Total 6"/8"

ADAPT PUMP™

Weight 11 lb.
Alternating Cycle Times 5/10/15 min.
Immersion Mode Standard
Microclimate Assist Standard

OTHER DETAILS

Alarms Power Failure/Low Pressure
Comfort Settings 8
Seamless Expandable Adapt Air EX™
Transport Cap. Standard

LEARN MORE

Ask your Agiliti representative or call 800-814-9389

