


# The Layers of NP Trio™

Three-Layer Advanced Foam Support Surface

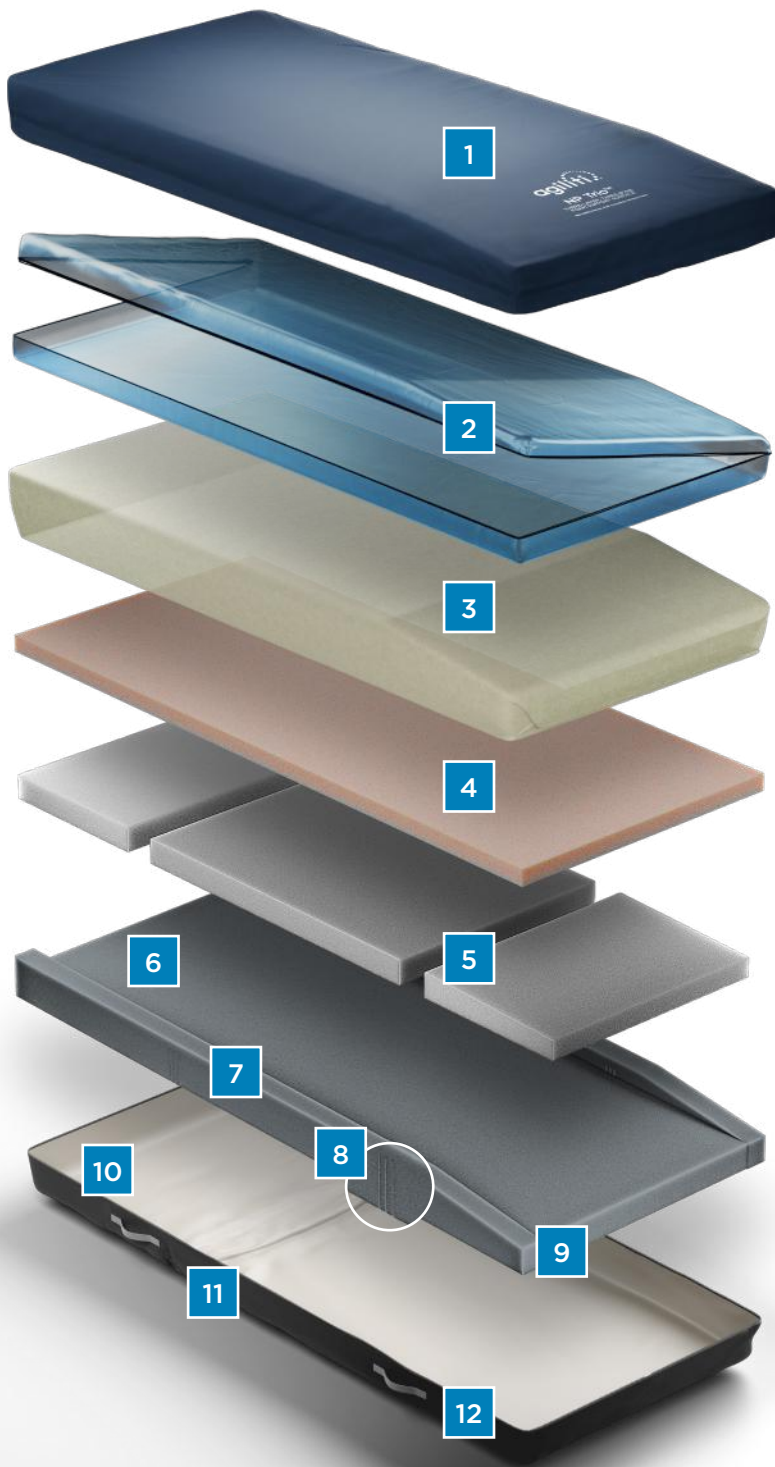

**Polycarbonate 4-way stretch top cover** (1) promotes immersion and envelopment. The breathable, chemically resistant performance fabric wicks away heat and moisture vapor and stands up to vigorous cleaning protocols.

**Core Shield™** (2) protects the internal components from fluid ingress and other damage, helping extend the life of the product.

The **halogen- and fiberglass-free fire barrier** (3) is safer for patients and caregivers—and better for the environment.

**CuPro™ copper- and gel-infused foam** (4) uses tiny particles of copper to help dissipate heat more efficiently—keeping the patient cool while reacting to temperature and weight.

A second layer of **Visco foam** (5) adds zoned support and pressure redistribution for the head, torso and heel sections.

A **third layer of foam** (6) provides extra support and a **firm foam edge** (7) adds stability during patient ingress/egress.

**Articulation grooves** (8) help the surface better contour to various bed angles, while the **sloped heel zone** (9) delivers additional off-loading capabilities.

The **bottom cover** (10) is RF-welded from a single panel of non-skid material to minimize the risk of fluid ingress and improve durability. Other details include an **expansion pleat** (11), which helps the surface contour to bed frame articulation and **color-coded handles** (12) to indicate year of manufacture.



**LEARN MORE**

Ask your Agiliti representative or call 800-814-9389