

# Adapt Air Pro™

Enhanced Immersion/Alternating Pressure Air Support Surface




8" Air Cells

Adapt Air Pro offers a combination of unique features designed for optimal pressure injury treatment and prevention. The 8" air cells deliver 30% more immersion compared to similar 6" support surfaces, while offering two clinically effective therapy modes: 1-in-2 alternating pressure or powered immersion. The support surface also features Microclimate Assist™ and AirSpace™ 3D mesh, which together provide constant airflow to the 4-way stretch top cover, helping keep patients cool and dry.



8" air cells provide 30% more immersion compared with 6" cells (Internal testing)

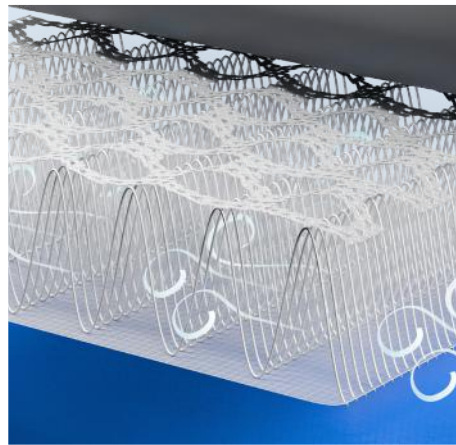

AirSpace 3D mesh maintains airflow while reducing friction and shear



Quiet, small-profile pump works on all Adapt Line™ support surfaces

## KEY SPECS (Standard/EX)

Max Patient Weight . . . . . 550 lb./1,000 lb.  
 Width . . . . . 35"/35"-48"  
 Length . . . . . 80"/82"-86"  
 Height: Air Cells/Total . . . . . 8"/8.5"

## ADAPT PUMP™

Weight . . . . . 11 lb.  
 Alternating Cycle Times . . . . . 5/10/15 min.  
 Immersion Mode . . . . . Standard  
 Microclimate Assist . . . . . Standard

## OTHER DETAILS

Alarms . . . . . Power Failure/Low Pressure  
 Comfort Settings . . . . . 8  
 Seamless Expandable . . . . . Adapt Air Pro EX™  
 Transport Cap. . . . . Standard

## LEARN MORE

Ask your Agiliti representative or call 800-814-9389

