

Adapt Air Pro™

Enhanced Immersion/Alternating Pressure Air Support Surface




8" Air Cells

Adapt Air Pro™ offers a combination of unique features designed for optimal pressure injury treatment and prevention. The 8" air cells deliver 30% more immersion compared to similar 6" support surfaces, while offering two clinically effective therapy modes: 1-in-2 alternating pressure or powered immersion. Targeted airflow microclimate management provides a constant flow of air to help keep patients cool and dry.



8" air cells provide 30% more immersion compared with 6" cells (Internal testing)

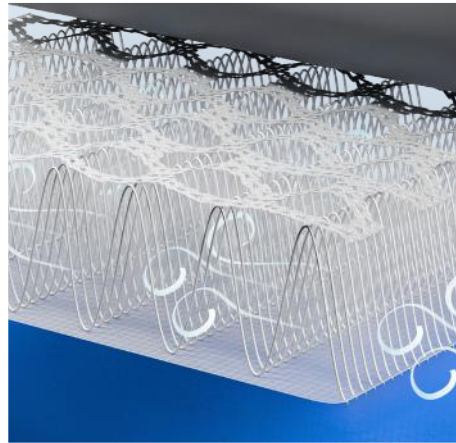

Targeted airflow microclimate management provides a constant flow of air to help keep patients cool and dry



Quiet, easy-to-use pump works on all Adapt Line™ support surfaces

KEY SPECS (Standard/EX)

Max Patient Weight 550 lb./1,000 lb.
 Width 35"/35"-48"
 Length 80"/82"-86"
 Height: Air Cells/Total 8"/8.5"

ADAPT PUMP™

Weight 11 lb.
 Alternating Cycle Times 5/10/15 min.
 Immersion Mode Standard
 Microclimate Management Standard


OTHER DETAILS

Alarms Power Failure/Low Pressure
 Comfort Settings 8
 Seamless Expandable Adapt Air Pro EX™
 Transport Cap. Standard

LEARN MORE

Ask your Agiliti representative or call 800-814-9389



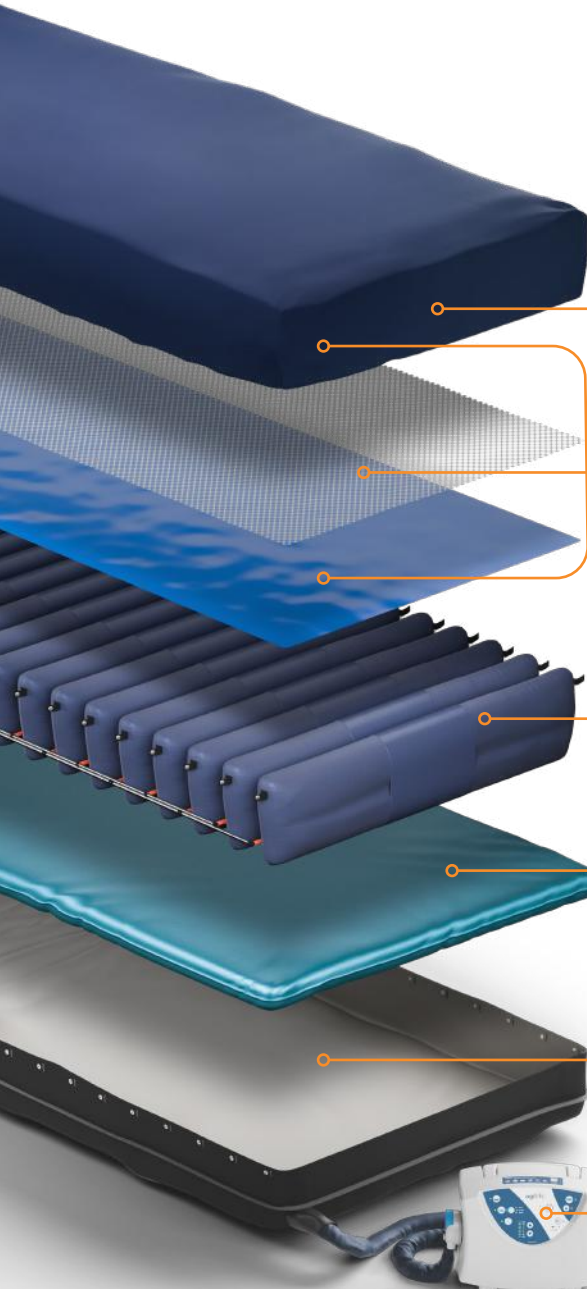

LAYER BY LAYER: ADAPT AIR PRO™

The **4-way stretch top cover** promotes immersion and envelopment; is breathable and chemically resistant, standing up to the toughest cleaning protocols.

Targeted airflow microclimate management combines the top cover fabric, Adapt Pump™ and **AirSpace™ mesh** to ensure non-stop airflow beneath the patient/surface interface — helping keep patients cool and dry.

Twenty 8-inch air cells provide alternating pressure or powered immersion. Unique integrated I-beams maintain vertical structure when weight is applied, leading to deeper immersion potential and lower shear forces.

Though the surface can maintain inflation without power, a **layer of foam** provides additional protection to help prevent bottoming out.

The **bottom cover** is RF-welded from a single panel of non-skid material to minimize the risk of fluid ingress and improve its durability.

The quiet, easy-to-use **Adapt Pump** is compatible across all Adapt Line™ support surfaces.

End-to-End Bed Services


Agility is the only manufacturer of therapeutic support surfaces with end-to-end services that support clinicians and patients, drive measurable outcomes and reduce cost.



Quick access to therapy options, both **rental and purchase**



Expedite time-to-therapy with **onsite management and logistics**

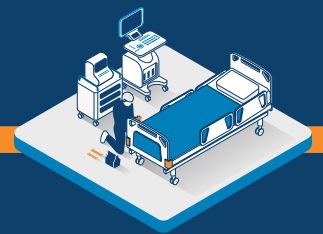

Optimize useful life with **comprehensive maintenance and repair**

Learn more at agilityhealth.com

PRD-0067G-0524 © 2024 Agility Health, Inc.

