

# Adapt Air Pro™

Enhanced Immersion/Alternating Pressure Air Support Surface




8" Air Cells

Adapt Air Pro™ offers a combination of unique features designed for optimal pressure injury treatment and prevention. The 8" air cells deliver 37% more immersion compared to similar 6" support surfaces, while offering two clinically effective therapy modes: 1-in-2 alternating pressure or powered immersion. Targeted airflow microclimate management provides a constant flow of air to help keep patients cool and dry.



8" air cells provide 37% more immersion compared with 6" cells (Internal testing)

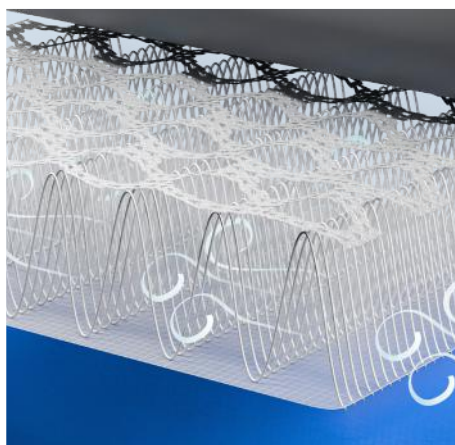

Targeted airflow microclimate management provides a constant flow of air to help keep patients cool and dry



Quiet, easy-to-use pump works on all Adapt Line™ support surfaces

## KEY SPECS (Standard/EX)

Max Patient Weight . . . . . 550 lb./1,000 lb.  
Width . . . . . 35"/35"-48"  
Length . . . . . 80"/82"-86"  
Height: Air Cells/Total . . . . . 8"/8.5"

## ADAPT PUMP™

Weight . . . . . 11 lb.  
Alternating Cycle Times . . . . . 5/10/15 min.  
Immersion Mode . . . . . Standard  
Microclimate Management . . . . . Standard


## OTHER DETAILS

Alarms . . . . . Power Failure/Low Pressure  
Comfort Settings . . . . . 8  
Seamless Expandable . . . . . Adapt Air Pro EX™  
Transport Cap. . . . . Standard

## LEARN MORE

Ask your Agiliti representative or call 800-814-9389



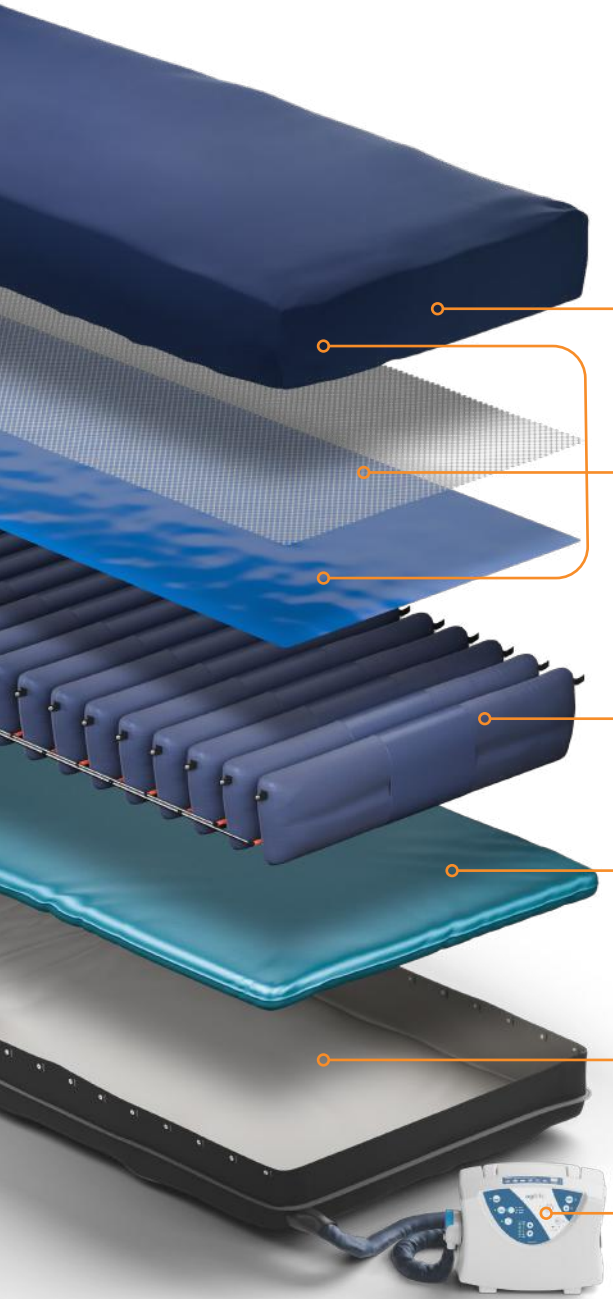

## LAYER BY LAYER: ADAPT AIR PRO™

The **4-way stretch top cover** promotes immersion and envelopment; is breathable and chemically resistant, standing up to the toughest cleaning protocols.

**Targeted airflow microclimate management** combines the top cover fabric, Adapt Pump™ and **AirSpace™ mesh** to ensure non-stop airflow beneath the patient/surface interface — helping keep patients cool and dry.

**Twenty 8-inch air cells** provide alternating pressure or powered immersion. Unique integrated I-beams maintain vertical structure when weight is applied, leading to deeper immersion potential and lower shear forces.

Though the surface can maintain inflation without power, a **layer of foam** provides additional protection to help prevent bottoming out.

The **bottom cover** is RF-welded from a single panel of non-skid material to minimize the risk of fluid ingress and improve its durability.

The quiet, easy-to-use **Adapt Pump** is compatible across all Adapt Line™ support surfaces.

## End-to-End Bed Services


Agiliti is the only manufacturer of therapeutic support surfaces with end-to-end services that support clinicians and patients, drive measurable outcomes and reduce cost.



Quick access to therapy options, both **rental and purchase**



Expedite time-to-therapy with **onsite management and logistics**

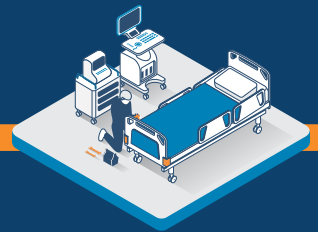

Optimize useful life with comprehensive **maintenance and repair**

Learn more at [agilitihealth.com](https://agilitihealth.com)

PRD-0067G-0924 © 2024 Agiliti Health, Inc.

