

Essentia™

Multi-Acuity Bed Frame
with SlideRail™



Essentia™ is a versatile, multi-acuity bed frame that can support a wide range of patients from Med-Surg to the ICU. Designed with high siderails to accommodate almost any therapeutic support surface, Essentia is a scalable platform that offers the flexibility for tailored care based on patient risk profile. New, innovative features – including the first-of-its-kind SlideRail™ – help lower caregiver risk and support patient mobility.



SlideRail: With up to 39" of unobstructed ingress/egress, SlideRail provides multiple secure handholds for patients getting in and out of bed.



Low Deck Height: In its lowest position, Essentia's deck is only 11.5" off the floor – among the lowest of any multi-acuity bed frame available.



RevTrend 20™: Start weight-bearing exercises sooner with a 20-degree reverse Trendelenburg position and reinforced footboard.

KEY SPECS

Safe Working Load 650 lb.
Max Patient Weight 550 lb.
Deck Low 11.5"
Deck High 30"
Deck Width 34"
Deck Length 80-88"
Head Angle 0° to 63°
Knee/Thigh Angle 0° to 39°
Foot Angle 0° to 9°
Trend/Reverse Trend -15° to 20°

OTHER DETAILS

Auxiliary Outlets Standard
Battery Backup Standard
Bed Exit Alarm Standard
CareMetrix™ Standard
Foot Extension (Electric) Standard
Hand Control Standard
MotionStop 360™ Standard
Nurse Call Connection Standard
PowerBrake™ w/AutoLock™ Standard
Protocol Timer Standard

SafeGlance™ Indicator Lights Standard
Scale System Standard
SlideRail™ Standard
Underbed Light Standard
USB Ports Standard


LEARN MORE

Ask your Agiliti representative or call 800-814-9389



Enhance Safety and Caregiver Workflow

Essentia™ simplifies how clinicians care for a broad range of at-risk patient populations. The bed's innovative features help protect caregivers and patients while enabling staff to more efficiently manage care.

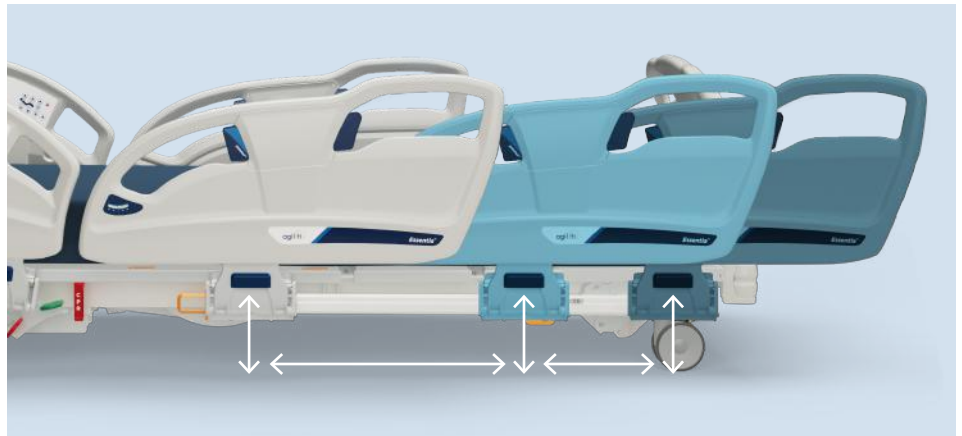

SlideRail™: Not only does it support patient safety during egress or ingress, SlideRail can also move out of the way — giving caregivers unobstructed access to provide bedside care.



SafeGlance™: Four indicator lights give busy caregivers a quick view of critical safety features, including bed exit alarm, protocol timer, bed low and brake status – features that save time and make it easier to manage the complex needs of a full patient unit.



PowerBrake™ with AutoLock™: With three locations on the bed, caregivers can easily set casters to brake, neutral or steer with the press of a button. And with AutoLock, the brakes lock after 60 seconds whenever the bed is plugged in.



Support Surface Options: A flat deck and higher siderails means Essentia is compatible with a wide range of support surfaces, including more than a dozen Agiliti-manufactured foam and air options.

End-to-End Bed Services


Agiliti is the only bed frame and therapeutic support surface manufacturer with end-to-end services that support clinicians and patients, drive measurable outcomes and reduce cost.



Quick access to product options, both **rental** and **purchase**



Expedite time-to-therapy with **onsite management and logistics**

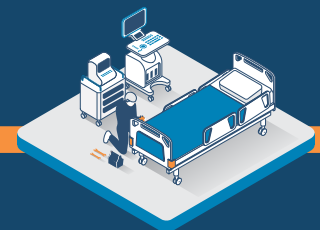

Optimize useful life with **comprehensive maintenance and repair**

Learn more at agilitihealth.com

PRD-0072B-0225 © 2025 Agiliti Health, Inc.

